

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

Convertible Chair Scale with BMI

MODEL MX308CHR

Unique convertible design gives you two scales in one: a wide body chair scale and a handrail scale. Conversion takes only seconds and requires no tools. First-class seating with 3" padded seat bottom and 25" span between armrests.

The "tilt & roll" design and robust wheels makes transportation of the scale effortless. Operation is simple and intuitive, even for the first time user.

The scale offers many **advanced features which come standard**, including Befour's proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an 'industry leading' 3-year warranty and are **Made in the USA**.



EMR READY

HANDS-FREE
WEIGHING

FEATURES & BENEFITS

Anti-Microbial Vinyl Covering: For sanitation and ease of cleaning

Slide-Out Footrest: For patient comfort

Extra Wide Chair: Accommodates most patient sizes

Integrated "Active" Handrail: When used as a stand-on scale, allows patient to grasp handrail for support while weight is accurately calculated

"Tilt & Roll" Design: Safely weigh patient when scale is stationary, but ready to move at any time

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

Convertible Chair Scale with BMI

MODEL MX308CHR

SCALE SPECIFICATIONS

SCALE HEIGHT 42"

CHAIR SCALE DIMENSIONS 29" W x 27" D x 42" H

SYSTEM WEIGHT 64 lbs

WEIGHT CAPACITY 500 lbs (220 kg)

RESOLUTION 0.1 lb (0.05 kg)

ACCURACY ± 0.1 lb (0.05 kg)

POWER SOURCE 4 D-cell alkaline batteries (incl)
*optional AC adapter available (P/N: 03049-08)

BATTERY LIFE 100,000 weighings

CONNECTIVITY USB serial or other options

DISPLAY READOUT 1" LCD digits

WARRANTY 3 years

SHIPPING SPECIFICATIONS


PACKAGED WEIGHT 81 lbs

PACKAGE DIM Box (45"x18"x10") & (30"x30"x25")

SHIPPING METHOD Freight

DIMENSIONS

Allow 1" perimeter space - scale should not touch wall

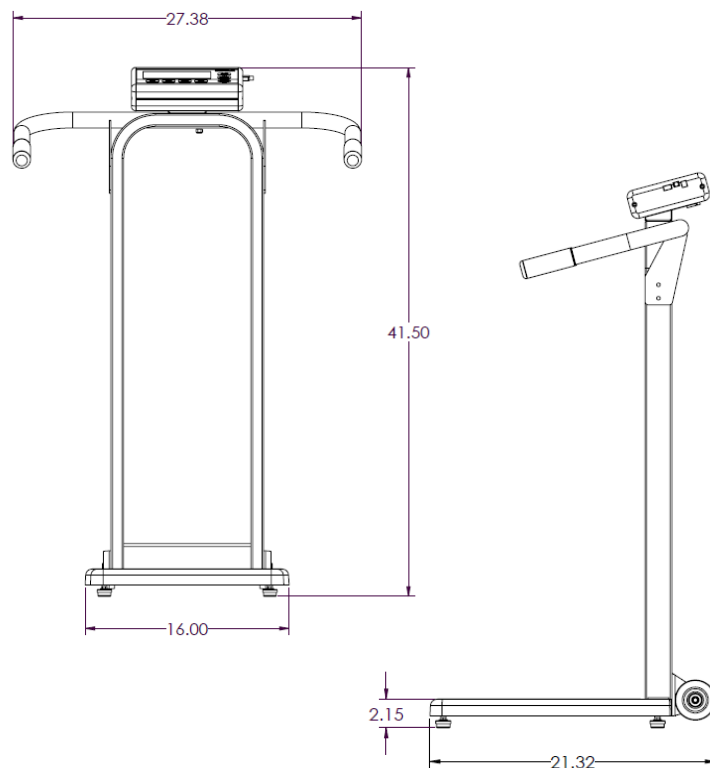

HANDRAIL SCALE CONVERSION



The MX308CHR Chair Scale quickly and easily converts to a HANDRAIL SCALE in seconds - all without the use of tools.

SLIDE-OUT FOOTREST

Slide-out footrest offers patient comfort while seated.



BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com