

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

Adjustable Width Portable Wheelchair Scale

MODEL MX420

The MX420 wheelchair scale is adjustable and portable, making it an excellent choice for the facility that occasionally needs to weigh patients in wheelchairs.

The dual rails adjust in width, from 16-3/4" to 35-1/2", to accommodate most wheelchairs. For mobile patients, this scale makes a **convenient chair scale** - simply place a chair on the rails and you are ready to weigh. This is popular option for visiting nurses, since it is lightweight and transports easily. The Display Console fits inside an **included carry case**.

Operation of the scale is simple and straight-forward, even for the first time user. The scale offers many **advanced features which come standard**, including Befour's proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an 'industry leading' 3-year warranty and are **Made in the USA**.



NOTE: The MX420 is recommended for use on hard flooring surfaces only.

FEATURES & BENEFITS

Adjustable-width rails: Conveniently slide to accommodate a variety of wheelchair widths

Ergonomic folding design: Folds in half for easy transport

High Capacity: 750 pounds

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

Adjustable Width Portable Wheelchair Scale

MODEL MX420

SCALE SPECIFICATIONS

CONSOLE MOUNTING Wall or tabletop

ADJUSTABLE DIMENSIONS 31" x 16 3/4" - 35 1/2"

SYSTEM WEIGHT 24 lbs

WEIGHT CAPACITY 750 lbs (340 kg)

RESOLUTION 0.1 lb (0.05 kg)

ACCURACY ± 0.1 lb (0.05 kg)

POWER SOURCE 4 D-cell alkaline batteries (incl)
*optional AC adapter available (P/N: 03049-08)

BATTERY LIFE 100,000 weighings

CONNECTIVITY USB serial or other options

DISPLAY READOUT 1" LCD digits

WARRANTY 3 years

SHIPPING SPECIFICATIONS


PACKAGED WEIGHT 31 lbs

PACKAGE DIMENSIONS Box (33" x 23" x 6")

SHIPPING METHOD UPS Package

DIMENSIONS

Allow 1" perimeter space - scale should not touch wall

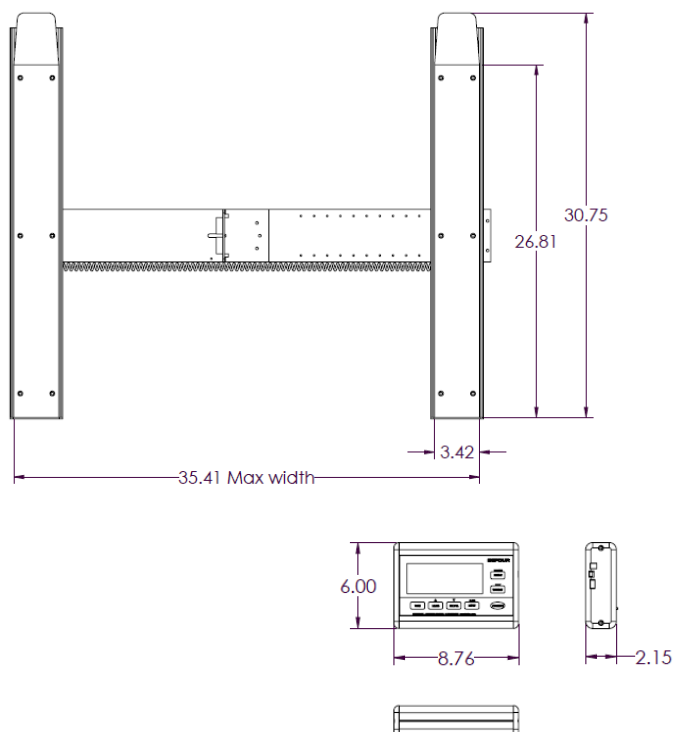

Versatile Console:

Mounting brackets (attached) allow positioning on wall or angled on tabletop



Scale folds up and display console stores in carry case (included) for easy transport and storage

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com