

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

Tilt & Roll Wheelchair Scale with Handrail

MODEL MX450 & MX450D

This versatile scale can be used as a wheelchair, handrail, or chair scale. Single ramp can be mounted on either side, while the dual ramp version allows for pass-through weighing. The ramps also fold up for use as a stand-on scale and the active handrail can be grasped during weighing.

The “tilt & roll” design and large wheels makes transportation of the scale effortless when moving from room to room.

The scale offers many **advanced features which come standard**, including Befour’s proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an ‘industry leading’ 3-year warranty and are **Made in the USA**.

FEATURES & BENEFITS:

Generous platform size: 28” x 32” platform accommodates most wheelchairs

Pivoting Console: Pivots 320° for optimal viewing

Integrated “Active” Handrail: When used as a stand-on scale, patients can grasp handrail for support and weight is accurately calculated

Simple Assembly: Minimal set-up required

“Tilt & Roll” Design: Safely weigh patient when scale is stationary, but ready to move at any time

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required


EMR READY



HANDS-FREE
WEIGHING



MX450

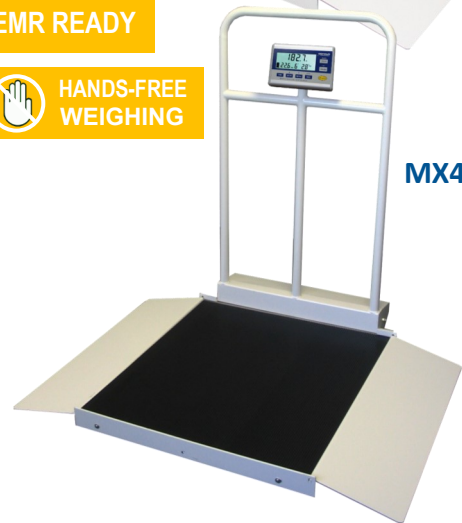

MX450D



DISPLAY CONSOLE

Display is used across all Befour’s medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, “Hands-Free” Weighing, etc.

‘ECO Weigh’ battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1” LCD digits

No AC outlet or recharging needed

Tilt & Roll Wheelchair Scale with Handrail

MODEL MX450 & MX450D

SCALE SPECIFICATIONS

SCALE HEIGHT 45"

PLATFORM DIMENSIONS 28"W x 32"D x 2"H

SYSTEM WEIGHT 65 lbs

WEIGHT CAPACITY 1000 lbs (450 kg)

RESOLUTION 0.1 lb (0.05 kg)

ACCURACY ± 0.1 lb (0.05 kg)

POWER SOURCE 4 D-cell alkaline batteries (incl)
*optional AC adapter available (P/N: 03049-08)

BATTERY LIFE 100,000 weighings

CONNECTIVITY USB serial or other options

DISPLAY READOUT 1" LCD digits

WARRANTY 3 years

SHIPPING SPECIFICATIONS




PACKAGED WEIGHT 78 lbs

PACKAGE DIMENSIONS Box (48" x 33" x 9")

SHIPPING METHOD Freight

DIMENSIONS

Allow 1" perimeter space - scale should not touch wall

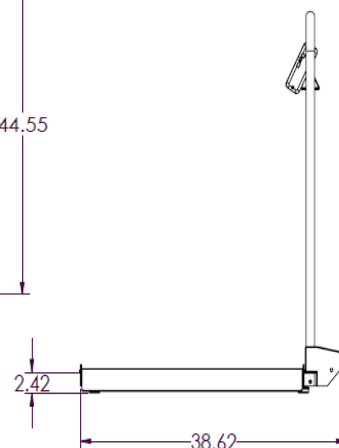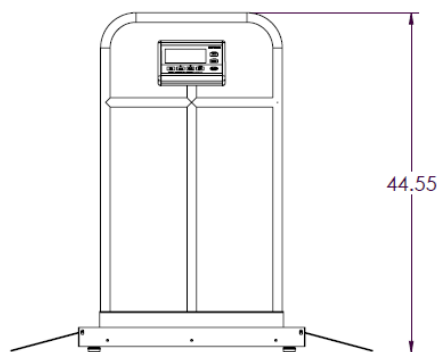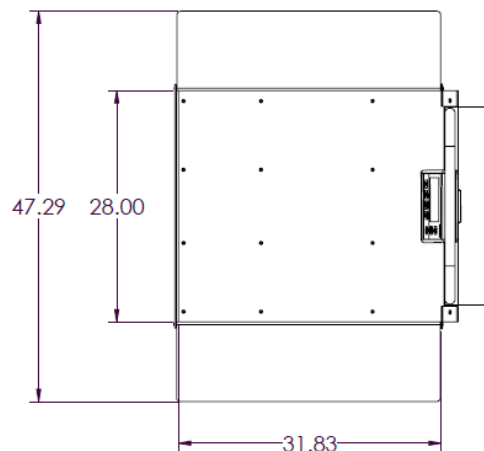

MX450 = Single Ramp
MX450D = Dual Ramp

TILT & ROLL

Convenient "Tilt & Roll" feature allows for easy transport to any location. Simply tilt scale back and gently roll forward.



OPTIONAL HEIGHT ROD (Mechanical)



The HTR450 Height Rod is not standard and needs to be ordered at time of scale as alterations to the scale need to be made during production

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com