

BEFOUR

PROFESSIONAL HEALTHCARE SCALES



MODEL MX905

The MX905 measurement station offers an attractive light-weight all-aluminum design, featuring a rugged cast platform, anodized column and Opti-Height Technology. This exclusive technology enables fast and accurate height measurements through state-of-the-art design and expert engineering. The 360° ergonomic hand grip and robust headpiece ensures the smoothest, simplest and safest height measurement process available in healthcare. The measurement slide encompasses Smart Position Technology meaning it 'stays put' where last left so height rod is always in the ready position for the next patient.

Operation is simple and intuitive, even for the first time user. This scale offers many **advanced features which come standard**, including Befour's proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an 'industry leading' **3-year warranty** and are **Made in the USA**.

EMR READY

HANDS-FREE WEIGHING

FEATURES & BENEFITS

Pin point accuracy of $\pm 0.1\text{lb}$

High Capacity: 750 pounds

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more

Optical Height Measurement: Simply lower the foldable headpiece to rest on top of head, press HEIGHT & measurement occurs instantaneously. Height can be viewed in selectable units and displays automatically

Built to Last: Rugged, non-corroding all-aluminum design weighs only 17 pounds.



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

Measurement Station Column Scale with BMI & Height Rod

MODEL MX905

SCALE SPECIFICATIONS

SCALE HEIGHT 42"

PLATFORM DIMENSIONS 16" W x 18" D x 1-3/4" H

SYSTEM WEIGHT 17 lbs

WEIGHT CAPACITY 750 lbs (340 kg)

RESOLUTION 0.1 lb (0.05 kg)

ACCURACY ± 0.1 lb (0.05 kg)

POWER SOURCE 4 D-cell alkaline batteries (incl)
*optional AC adapter available (P/N:03049-08)

BATTERY LIFE 100,000 weighings

CONNECTIVITY USB serial or other options

DISPLAY READOUT 1" LCD digits

HEIGHT ROD RANGE 25-94" / (~63-238 cm)

HEIGHT ROD ACCURACY ± 0.1 in (0.2 cm)

HEIGHT DISPLAY FORMATS In, Ft & In, Cm

WARRANTY 3 years

SHIPPING SPECIFICATIONS

PACKAGED WEIGHT 29 lbs

PACKAGE DIMENSIONS Box (45" x 18" x 10")

SHIPPING METHOD UPS Package


DIMENSIONS

Allow 1" perimeter space - scale should not touch wall



OPTI-HEIGHT DIGITAL HEIGHT ROD

Height rod is both digital and integrated into scale for automatic BMI calculation. Unique design has extended measuring range with no exposed cables and rotates 360° for easy positioning of headpiece.

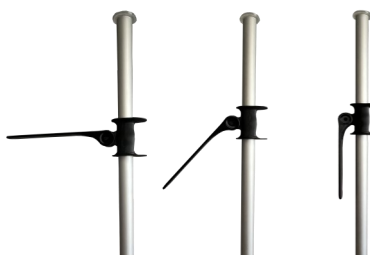

Measurement Slide has **Smart Position Technology** and stays where last left so it's always ready to measure.

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com