

BEFOUR

PROFESSIONAL HEALTHCARE SCALES

THE Exam Room Scale

MODEL MX877

Befour's **PREMIER** physician scale - THE Exam Room Scale - offers every feature you'd want in a stationary patient scale.

Stability, safety and technology are at the forefront of engineering. An innovative 'dual axis tip-proof' design and integrated "active" handrail - which comes standard - provide unmatched safety. This free-standing, rock-solid, scale incorporates all-metal construction to deliver robust quality in a sleek, space-saving package.

The model MX877 of THE Exam Room Scale features a 750 pound capacity, Befour's proprietary MotionLock technology, and ECO-Weigh power management system providing more than 100,000 weight readings on a single set of batteries. The MX877 includes a 360° rotatable integrated digital height rod mounted on the right side of the scale. Should your room layout dictate, the scale can be specified with the digital height rod located on the left side of the display console (MX877L).

All Befour scales are **built to last**, backed by an 'industry leading' **3-year warranty** and are **Made in the USA**.

*Model **MX877** pictured with standard height rod positioning

FEATURES & BENEFITS

Integrated "Active" Handrail: A **STANDARD** feature, allowing patients to grasp handrail for support while weight is accurately calculated

'Dual Axis Tip-Proof' Design: Unique, free-standing design eliminates all front-to-back and side-to-side tipping for unmatched stability and patient security

Display Console: LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more

Digital Height Measurement: When ordered as an MX877, simply press HEIGHT & raise the foldable headpiece. Height can be viewed in selectable units & displays automatically



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required

EMR READY



HANDS-FREE
WEIGHING



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

THE Exam Room Scale

MODEL MX877

SCALE SPECIFICATIONS

SCALE FOOTPRINT 16" W x 22" D x 46" H

PLATFORM DIMENSIONS 16" W x 18" D x 2" H

SYSTEM WEIGHT 41-44 lbs

WEIGHT CAPACITY 750 lbs (340 kg)

RESOLUTION 0.1 lb (0.05 kg)

ACCURACY ± 0.1 lb (0.05 kg)

POWER SOURCE 4 D-cell alkaline batteries (incl)
*optional AC adapter available

BATTERY LIFE 100,000 weighings

CONNECTIVITY USB serial or other options

DISPLAY READOUT 1" LCD digits

HEIGHT ROD RANGE 38-88" / (97-225cm)

HEIGHT ROD ACCURACY ± 0.1 in (0.2 cm)

HEIGHT DISPLAY FORMATS In, Ft & In, Cm

WARRANTY 3 years



SHIPPING SPECIFICATIONS

PACKAGED WEIGHT ~ 45-48 lbs

PACKAGE DIMENSIONS ~ Box (48 x 20" x 10")

SHIPPING METHOD UPS Package

DIMENSIONS

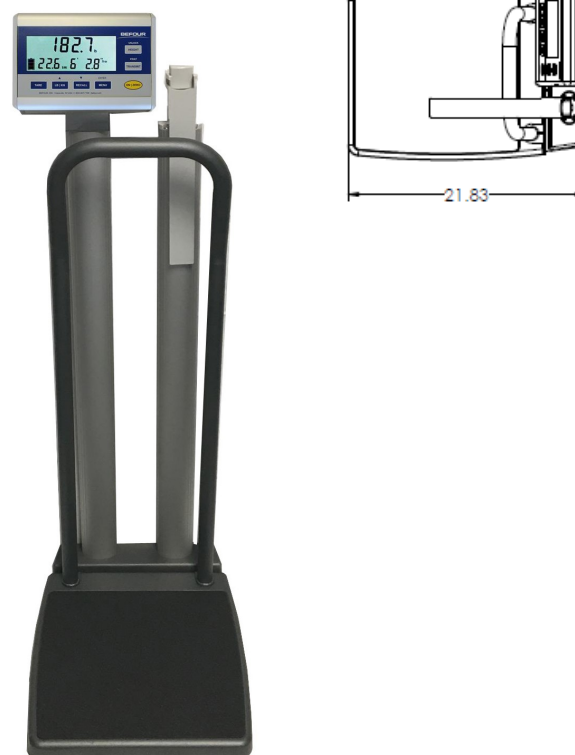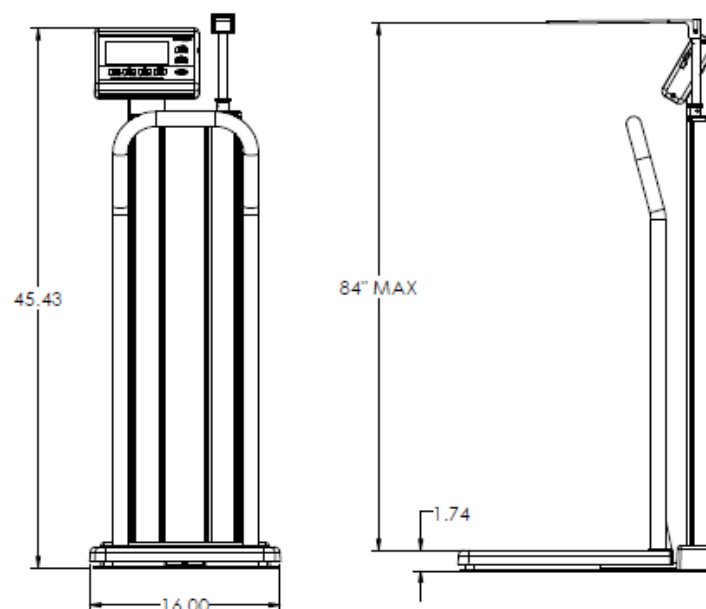

DUAL AXIS TIP-PROOF DESIGN *

Befour's exclusive 'dual axis tip-proof' engineering offers rock-solid stability, in a free-standing design, producing an incredibly secure scale.

* Design reduces chance of tipping under normal usage.



BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com