

BEFOUR


SPORTS & FITNESS SCALES

Portable Battery-Powered Bariatric Scale with BMI

Model PS-7700

The PS-7700 portable scale is ideal for bariatric professionals, home health care, athletic trainers, weight-loss clinics, and weight-lifting competitions due to its generous 24" x 20" platform, 900 lb weight capacity, and built-in body mass index calculator.

The low-profile platform sits just 2" off the ground, making it easy to step on for individuals of all sizes. Battery powered scale-no outlet or re-charging needed!

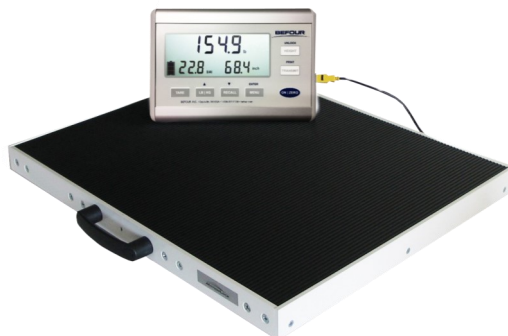

FEATURES & BENEFITS:

Body Mass Index: Dual-mode BMI calculator allows user-or assistant-operation and provides a quick and easy BMI calculation, programmable clothing allowance built-in

Pin point accuracy of ± 0.1 lb.

High Capacity: 900 lbs.

Fully assembled: Arrives ready to use; simply insert the batteries (included)

Display Console Set-up: Quick insert/detach connections

Push Button Calibration: Set calibration via keypad entry without the need of tools

Rugged and durable all-metal construction

Lightweight and Easy to Carry: Build-in handle and recessed storage cavity securely holds the display console during storage and extended transport

MOTIONLOCK Accurately calculates the weight of an unsteady user

Plug & Weigh Exclusive cable system can be easily replaced if damaged

SPECIFICATIONS:

Platform Dimensions.....	20" x 24" x 2"
System Weight	25 lbs
Weight Capacity.....	900 lbs (225 kg)
Resolution.....	0.1lb (0.05 kg)
Accuracy.....	± 0.1 lb (0.05kg)
Power Source.....	D C-cell alkaline batteries
Battery Life.....	100,000+ weighings
Display Readout.....	1" LCD digits
Warranty.....	3 years

DISPLAY CONSOLE FEATURES



- **Displays weight, height and BMI** on a single screen
- **Programmable User Settings:** can be set to auto-off (12 seconds to 25 minutes), always-on, etc.
- **LB/KG Lockout:** Change setup to display only your desired unit of measure
- 3-segment battery life indicator
- Large 1" LCD digits
- No AC outlet or recharging needed (batteries last 100,000+ weighings)
- Optional AC Adapter available: P/N 03049-06
- Auto zero, recall, and reweigh

Built-in handle and recessed storage cavity securely holds the display console during extended transport.



Befour, Inc.

102 Progress Drive · Saukville, WI 53080 · P: 800-367-7109 · F: 262-284-5966
Befour.com · sales@befour.com