

BEFOUR

SPORTS & FITNESS SCALES

Model WH-1062

Befour's WH-1062 weight and height measurement scale features our Opti-Height Optical Stadiometer System. The WH1062 quickly determines weight, height and BMI. With a weight capacity of 650 lbs., +/- 0.1lb accuracy, and MotionLock technology, which provides an accurate weight of unsteady users, this scale is an exceptional value.

The incredible battery life provides years of use on a single set of batteries—no outlet or recharging needed.

FEATURES & BENEFITS:

Pin point accuracy of ± 0.1 lb.

High capacity: 650 lbs.

Body Mass Index: Automatically calculates and displays BMI

Push Button Calibration: Set calibration via keypad entry without the need of tools

Rugged and durable all-metal construction

Display Console: LB/KG lockout, programmable settings, BMI output, and more highlighted in right-hand column

Optical Height Measurement: Simply lower the foldable headpiece to rest on top of head, press HEIGHT & measurement occurs instantaneously. Height can be viewed in selectable units and displays automatically

 Accurately calculates the weight of an unsteady user

SPECIFICATIONS:


Height Rod Range.....	25-94" (63-238cm)
Platform Dimensions.....	16" x 18" x 1 3/4"
System Weight.....	17 lbs
Resolutions.....	0.1 lb (0.05 kg)
Accuracy.....	± 0.1 lb (0.05 kg)
Power Source.....	4 D-cell Alkaline Batteries
Battery Life.....	100,000 Weighings
Display Readout.....	1" LCD digits
Warranty.....	3 years



DISPLAY CONSOLE FEATURES



- **Displays weight, height and BMI** on a single screen
- **Programmable User Settings:** can be set to auto-off (12 seconds to 25 minutes), always-on, etc.
- **LB/KG Lockout:** Change setup to display only your desired unit of measure
- 3-segment battery life indicator
- Large 1" LCD digits
- No AC outlet or recharging needed (batteries last 100,000 weighings)
- Optional AC Adapter available: PN 03049-06
- Auto zero, recall, and reweigh

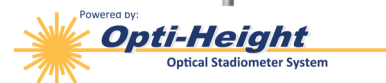

Height rod is both digital and integrated into scale for automatic BMI calculation. Unique design has extended measuring range with no exposed cables and rotates 360° for easy positioning of headpiece.



Befour, Inc.

102 Progress Drive · Saukville, WI 53080 · P: 800-367-7109 · F: 262-284-5966
Befour.com · sales@befour.com