

# BEFOUR

PROFESSIONAL HEALTHCARE SCALES



Powered by:

## Opti-Height

Optical Stadiometer System

### ADA COMPLIANT MODEL MX908

Compact, lightweight, and portable, this measuring station is designed to maximize your healthcare investment and is engineered with Opti-Height Technology. This exclusive feature enables fast and accurate height measurements through state-of-the-art design and expert engineering. The 360° ergonomic hand grip and robust headpiece ensures the smoothest, simplest and safest height measurement process available in healthcare. The measurement slide encompasses *Smart Position Technology* meaning it 'stays put' where last left so height rod is always in the ready position for the next patient.

Transporting this scale is effortless due to the *"tilt & roll" design and robust wheels* and operation is simple and intuitive.

The scale offers many *advanced features which come standard*, including Befour's proprietary MotionLock technology and ECO Weigh battery system for years of use and allows the scale to be used anywhere.

All Befour scales are **built to last**, constructed of anodized aluminum, backed by an 'industry leading' **3-year warranty** and **Made in the USA**.



## FEATURES & BENEFITS

**Integrated "Active" Handrail:** Allows patient to grasp handrail for support while weight is accurately calculated

**Generously Sized Wrap-Around Handrail:** Full 25" spacing between handgrips. Unique, heavy-duty, design does not interfere with foot positioning and adds security and stability for patients of all sizes

**"Tilt & Roll" Design:** Safely weigh patient when scale is stationary, but ready to move at any time

**Display Console:** LB/KG lockout, programmable settings, BMI output, auto zero, recall, reweigh ability, and more

**Optical Height Measurement:** Simply lower the foldable headpiece to rest on top of head, press HEIGHT & measurement occurs instantaneously. Height can be viewed in selectable units and displays automatically



MotionLock automatically facilitates quick & accurate weight, even for the most unsteady user – no intervention of the care provider or patient required



DISPLAY CONSOLE

**Display is used across all Befour's medical scales:** Seamlessly move from one model to another with ease

**Displays weight, height, and BMI** on a single screen

**EMR connectivity** via serial output or other options

**Programmable User Settings** can be set to auto-off, "Hands-Free" Weighing, etc.

**'ECO Weigh' battery system** provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

**LB/KG Lockout:** Toggle between or set to display only your desired unit of measure

3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

# Measurement Station Tilt & Roll Handrail Scale

## MODEL MX908

### SCALE SPECIFICATIONS


|                        |                                    |
|------------------------|------------------------------------|
| SCALE HEIGHT           | 42"                                |
| PLATFORM DIMENSIONS    | 16" W x 18" D x 1-3/4" H           |
| SYSTEM WEIGHT          | 33 lbs                             |
| WEIGHT CAPACITY        | 750 lbs (340 kg)                   |
| RESOLUTION             | 0.1 lb (0.05 kg)                   |
| ACCURACY               | ± 0.1 lb (0.05 kg)                 |
| POWER SOURCE           | 4 D-cell alkaline batteries (incl) |
| BATTERY LIFE           | 100,000 weighings                  |
| CONNECTIVITY           | USB serial or other options        |
| DISPLAY READOUT        | 1" LCD digits                      |
| HEIGHT ROD RANGE       | 25-94" / (~63-238 cm)              |
| HEIGHT ROD ACCURACY    | ± 0.1 in (0.2 cm)                  |
| HEIGHT DISPLAY FORMATS | In, Ft & In, Cm                    |
| WARRANTY               | 3 years                            |

### SHIPPING SPECIFICATIONS

|                    |                       |
|--------------------|-----------------------|
| PACKAGED WEIGHT    | 44 lbs                |
| PACKAGE DIMENSIONS | Box (45" x 18" x 10") |
| SHIPPING METHOD    | UPS Package           |


### DIMENSIONS

Allow 1" perimeter space - scale should not touch wall

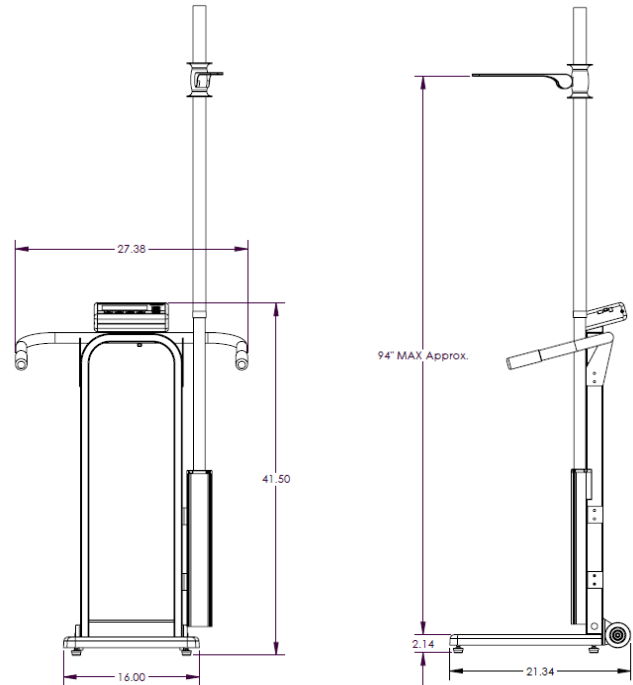


### OPTI-HEIGHT DIGITAL HEIGHT ROD

Height rod is both digital and integrated into scale for automatic BMI calculation. Unique design has extended measuring range with no exposed cables and rotates 360° for easy positioning of headpiece.

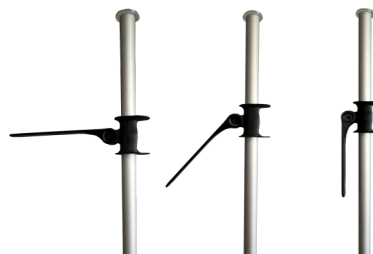

Measurement Slide has *\*Smart Position Technology\** and stays where last left so it's always ready to measure.

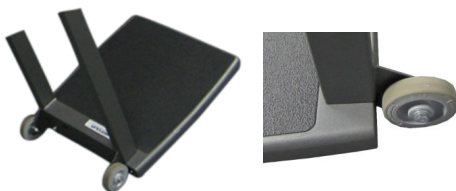

### TRANSPORT WHEELS

To move to a new location, simply tilt the scale back on its robust 3" wheels and roll away. The scale fits through all standard doorways.

**BEFOUR**

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com