

Opti-PED

Pediatric Exam Device

A complete pediatric exam station that safely measures weight and length of neonates to toddlers...in seconds!

SCALE + INFANTOMETER + EXAM SURFACE

- Weigh active infants quickly and precisely - 10 gram accuracy
- Displays digital length at the push of a button - 1mm accuracy
- Available as a stand-alone unit or with optional table and storage



DISPLAY CONSOLE

Display is used across all Befour's medical scales: Seamlessly move from one model to another with ease

Displays weight, height, and BMI on a single screen

EMR connectivity via serial output or other options

Programmable User Settings can be set to auto-off, "Hands-Free" Weighing, etc.

'ECO Weigh' battery system provides 100,000 weighings per set of (4) D-cell alkaline batteries (included)

LB/KG Lockout: Toggle between or set to display only your desired unit of measure


3-segment battery life indicator

Large 1" LCD digits

No AC outlet or recharging needed

Opti-PED

Pediatric Exam Device

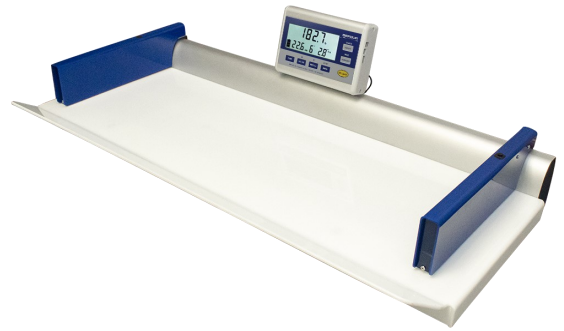

Opti-PED delivers unsurpassed safety, accuracy, capacity, speed, and ease of cleaning. With its “no-touch” display console, Opti-PED permits health professionals to focus attention where it should be - on the infant. Operation is simple. Place infant on acrylic tray and weight instantly displays. To measure length, there is no need to reposition the infant. Effortlessly guide the dual paddles into position, press paddle button and length digitally appears. All measurement values stay displayed even after the infant is removed.

FEATURES / SPECIFICATIONS

Self-Contained: Arrives fully assembled and ready to use.

Neonates to Toddlers: With 10 gram accuracy & 27 kilogram (60 pound) capacity, Opti-PED can handle all your pediatric exam needs.

One-Position Measurement: Infant remains in same position for the entire weight and length measurement sequence. Tray has gently sloping front for easy access to infant for examination.

LB/GM Lockout: Quickly change setup to only weigh in selected unit of measure: kg, lb, or lb/oz (weight) and decimal cm, decimal inch or fractional inch (length).



Accurately calculates the weight of active babies.



Model TBL282 - Optional Table

Opti-PED becomes a portable pediatric exam system with the TBL282

- Easy – clean stainless steel construction, 24” x 48” top surface
- Multiple storage drawers
- 500 pound capacity
- Provides ergonomic 35” to 36” working height for Opti-PED tray
- Includes 14” exam paper holder & cutter bar
- Includes swivel casters (2 locking) for maximum security



PRODUCT SPECIFICATIONS

Exam Tray Dimensions.....	39”w x 18”d x 12”h
Tray Working Height.....	3” above countertop surface
System Weight.....	37 lbs (17kg)
Scale Capacity.....	60 lbs (27 kg)
Scale Resolution.....	0.01 lb / 0.01 oz (0.01kg)
Scale Accuracy.....	± 0.02 lb / 0.3 oz (0.01 kg)
Infantometer Range.....	9 to 39 in (23 to 97 cm)
Infantometer Resolution.....	0.1 in / 1/8 in (0.1cm)
Infantometer Accuracy.....	± 0.1 in / 1/8 in (0.1cm)
Power Source	4 D-cell alkaline batteries
	*optional AC adapter available (P/N:03049 –08)
Battery Life (Always On).....	Approximately 1 year
Display.....	1” LCD
Warranty.....	3 years

SHIPPING SPECIFICATIONS

Packaged Weight.....	45lbs
Package Dimensions.....	48”L x 27 ”W x 14”H

BEFOUR

Befour Inc. | 102 Progress Drive, Saukville, WI 53080 | 800.367.7109 | f: 262.284.5966 | befour.com | sales@befour.com